




## Product description:

The program is intended to collect both partial weighing procedures and packed product control complete reports from Radwag checkweighers.

The application can serve many balances at the same time, enabling automatic segregating and saving loaded weighing procedures in separate databases. It also enables reports generating and exporting data to text format.

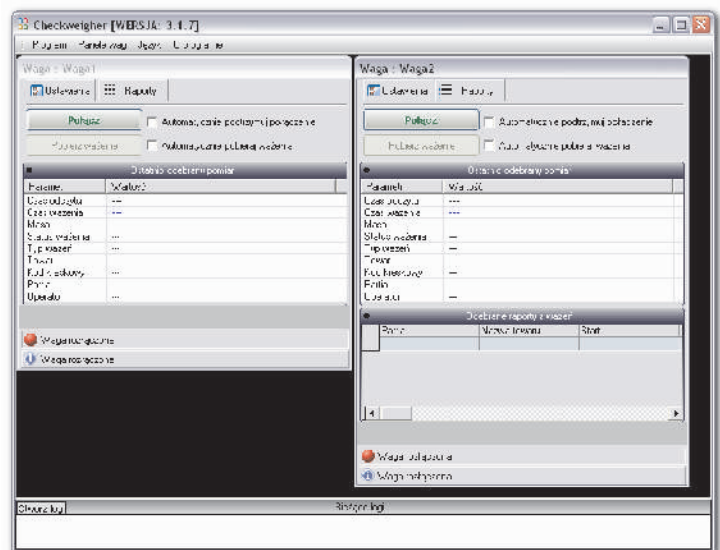

## Designed for scales:

- dynamic/control, based on TMC head

## Functions:

- servicing many balances at the same time
- program operating in automatic or manual mode, measurements collecting from balance heads and sending them to local databases
- collected measurements segregating inside database (MS Access), giving consideration to time and measurement category
- measurements (databases) exporting to CSV text format, with semicolon as records separator
- communication servicing through port RS-232 and Ethernet (TCP)

## Computer requirements:

- processor 1,2 GHz
- free 500 MB on hard drive
- RAM 256 MB (recommened 512 MB)
- operation system Windows 2000/XP

## Language versions:

Polish, English, Russian.

\* NOTE! RADWAG does not sell software "in box". Boxes presented on our Web-site are of virtual characteristic only, and they are not accessible with purchased software. The customer receives purchased software on a CD in special protective cover and the user manual as a brochure.