

# JEWELLERY BALANCES AS




Analytical balances AS series are offered with following capacities: 140 g, 220 g and 320 g with readability of 0,2 mg.

They are equipped with anti breeze draught shield with sliding glass doors and backlit LCD display.

These balances have system of internal automatic calibration (C2 version) activated on time and temperature changes. Balances are equipped with possibility of weighing loads outside the main weighing platform (so called under hook weighing). This means of mass measuring is an alternative for loads with non-standard dimensions and shapes and those which create magnetic field. Under hook weighing is also applied for density determination procedures.

- Counting pieces
- Density determination
- Percentage
- Summing function

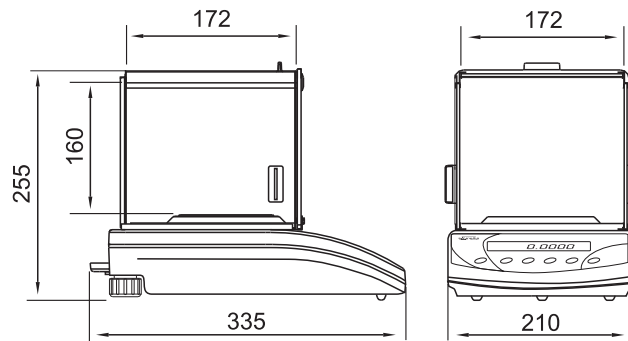

## Technical data:

	AS 140/C/1/CT AS 140/C/2/CT *	AS 220/C/1/CT AS 220/C/2/CT *	AS 320/C/1/CT AS 320/C/2/CT
Max capacity	140 g / 700 ct	220 g / 1100 ct	320 g / 1600 ct
Minima load	0,004 g / 0,02 ct	0,004 g / 0,02 ct	0,004 g / 0,02 ct
Readability	0,0002 g / 0,001 ct	0,0002 g / 0,001 ct	0,0002 g / 0,001 ct
Tare range	-140 g / -700 ct	-220 g / -1100 ct	-320 g / -1600 ct
Working temperature	+15° - +30° C		
Reapatability	0,0002 g / 0,001 ct	0,0002 g / 0,001 ct	0,0004 g / 0,002 ct
Linearity	±0,0002 g / 0,001 ct	±0,0002 g / 0,001 ct	±0,0004 g / 0,002 ct
Sensitivity drift	2 ppm/°C in temperature +15° - +35°C		
Power supply	230V / 11V AC or 120V / 11V AC		
Calibration	internal (C2 version) or external (C1 version)		
Pan size	ø 85 mm		

\* verified balances

## Additional equipment:

ANTI-VIBRATION TABLE	SUITCASE FOR BALANCE
"KAFKA" THERMAL PRINTER	LCD ADDITIONAL DISPLAY
NEEDLE PRINTER "EPSON"	EXTERNAL RECHARGEABLE BATTERY PACK
DENSITY DETERMINATION KIT FOR SOLIDS AND LIQUIDS	COMPUTER SOFTWARE: PW-WIN 2005, RAD-KEY 2000, REC-FS 2000
RACK FOR UNDER HOOK WEIGHING	

## RADWAG USA L.L.C.

19599 NE 10th Ave., Bay G, • North Miami Beach, FL, 33179 • USA • Tel: 1-305-651-3522 • Fax: 1-305-651-3523 • e-mail: office@radwagusa.com • website: www.radwagusa.com